




MELTING POINT APP AC-361



MELTING POINT AC-361

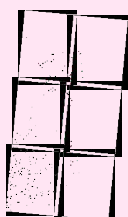

COOLING PLUG



DIGITAL MODEL AC-361B

SPECIFICATIONS

ITEM No. 1 MELTING POINT APPARATUS (CAT. No. AC-361)

- Body : All Aluminum Body
- Usage : For Determination of Single or Mized Melting Points.
- Temperature : Room Temp. to 350°C
- Heater Block : All Block accepts 3 capillary tubes & Mercury thermometer
- Heating : The above Block is heated by two elements clumped to the sides
- Illumination & Viewing : The built-in lamp view through a bocular lens
- Temp. Control : A Rotary Thermostat

- Optional: « Cooling plug for water circulation.
- : Thermometers 300°C

ITEM No. 2 MELTING POINT APPARATUS DIGITAL (CAT. No. AC-361B) digi

- Features : All same as above but having Built in Probe and Digital display for temperature attained while watching solid turning to liquid

AC-361 Capillary thermostat

AC-361 Digital PIC Controlled Model

As above but digital temperature controller cum indicator
 Range Temperature 5° above room to 100°C Least count)
 0.1°C with SSR Relay

AC-361 is Our new Catalogue Model of AH-361



ABRON EXPORTS

SCIENTIFIC INSTRUMENTS COMPANY

We Proud to be AN ISO 9001:2008 Company

